

NEW
TRACHWAY[®]
VIDEO INTUBATION SERIES

Trachway Video Intubation Series offers
total solution for clinical tracheal intubation



Video Light Stylet

TRACHWAY[®]
VIDEO INTUBATION SERIES





Features

- Able to work inside oral space without stimulation
- Bendable and not easily deformed
- Visualization nasal Intubation
- Trachlight function assist

Operation board

Recording function

Snapshot 

Video record 




Power source

Click-in connection
Quick switch among blades,
stylets, and flexible stylet



Malleability

Bend for difficult airway case
and nasal intubation

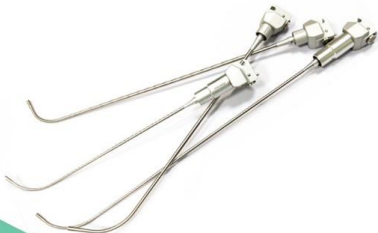

Advantage

- Patented guiding red LED
- 304 stainless steel
- Submersible for disinfection
- Seamless seal and foggy free



Quick switch connector

Click-in design for quick switch
among blade, style,
and soft style systems



Four Sizes

TVI-4102 Ø 5.0mm L: 345mm
-General ETT (5.5 above)
-Reinforce ETT

TVI-4105 Ø 5.0mm L: 420mm
-35 Fr. DLT, Preform ETT

TVI-4106 Ø 4.0mm L: 300mm
-Pediatric ETT (4.5 above)

TVI-4108 Ø 4.0mm L: 420mm
-32 Fr. DLT, Preform ETT



Video Flexible Stylet

TRACHWAY[®]
VIDEO INTUBATION SERIES

Features

Visualization tube exchanger

Double lumen ETT placement and confirmation

Percutaneous tracheostomy procedure monitoring



Flexible Tube

Not easy to break



Plug in & out connector

Fast and efficient when implement tube exchange



Two Sizes

TVI-4301 OD 5.0mm Lenth 750mm

TVI-4302 OD 4.0mm Lenth 750mm



Video Laryngoscope



Features

- Conventional blade design
- No learning curve
- Improve one Cormack & Lehane grade
- Submersible blades



Release latch

Separation monitor and handle for high chest patient and teaching purpose



Real time video display
Wireless transmission for real time video



Durability

- Submersible for disinfection
- Seamless seal
- Foggy free



Advantage

- Conventional blade design
- Same skill with non direct view
- 316 stainless steel



Five Sizes

- TVL-4220 Mill#0
-W: 11.5mm L: 75mm
- TVL-4221 Mill#1
-W: 11.5mm L: 102mm
- TVL-4232 MAC#2
-W: 12mm L: 100mm
- TVL-4233 MAC#3
-W: 15mm L: 130mm
- TVL-4234 MAC#4
-W: 15mm L: 135mm



Quick switch connector

Click-in design for quick switch among blade, style, and soft style systems



Specification



▶ Trachway® Handle Set TVI-4000



Adjustable angle
Thumb operation
60° upward
45° downward



- Video output port**
- Micro SD card slot**
- Mini USB Charging port**



▶ TVI-4051 3.5" Monitor

Display	3.5 " TFT
Power Source	3.7V Li-ion rechargeable battery
Charger	AC adaptor with DC5V, 2.4A output
Casing material	ABS/PC
Weight	174g
Dimensions	W 90 x H 110 x D 88 (mm)



▶ TVI-4050 Rechargeable Handle

Power consumption	3.7V Li-ion rechargeable battery
Full Charge Operation	≥ 4hrs
Luminance adjustment	5 Levels
Casing material	ABS/PC
Weight	115g
Dimensions	W 32 x H 111 x D 39 (mm)



▶ Video Light Stylet

Light Source	Surround LEDs Guiding red LED (for Ø 5.0)
Material	S.S. 304
Focus Range	5-50mm
Field of view	FOV (D): 83°

- TVI-4102 Ø 5.0 L 345mm
- TVI-4105 Ø 5.0 L 420mm
- TVI-4106 Ø 4.0 L 300mm
- TVI-4108 Ø 4.0 L 420mm



▶ Video Flexible Stylet

Light Source	Surround LEDs
Material	Medical Grand PU
Focus Range	5-50mm
Field of view	FOV (D): 83°

- TVI-4301 Ø 5.0 L 650mm
- TVI-4302 Ø 4.0 L 750mm



▶ Video Flexible Blade

Light Source	LEDs
Material	S.S. 316
Focus Range	5-50mm
Field of view	FOV (D): 83°

- TVL-4220 MILL #0 Neonatal
- TVL-4221 MILL #1 Infant
- TVL-4232 MAC #2 Pediatric
- TVL-4233 MAC #3 Adult
- TVL-4234 MAC #4 Adult

Accessorize



▶ **TVI-2060** Tube resistor #5.0
(for video light stylet OD 5.0)



▶ **TVI-2061** Tube resistor #4.0
(for video stylet OD 4.0)



▶ **TVI-4056-B** Aluminum carrying case
L 508 W 270 H 100 (mm)
Weight 1.5kg

MANUFACTURER



榮晟醫療器材股份有限公司
Markstein Sichte Medical Corp.

2F., No.91, Ganhe Rd.,
Xitun Dist., 40751
Taichung,
Taiwan

T: 886-4-2315-3088
F: 886-4-2315-2589
E: fd@msst.com.tw
www.msst.com.tw



DISTRIBUTOR

